



8 Well chambered cover Glass with #1.5 high performance cover glass - 57mm x 25mm base

- Technical specs
- Dimension diagram
- Products with coating
- Cited publications



8 Well Chambered Cover Glass with #1.5 high performance cover glass (0.170±0.005mm), with lid, sterilized. Designed for high resolution imaging such as confocal microscopy.

Coverslip :
#1.5H » [view coverslip specs](#)

Catalog # :
C8-1.5H-N , [request a free sample](#) or [Get a quote](#)

Packing :
48/case

Price :
\$274.00 USD/case

1 case ▼

[+ Add to Cart](#)

Availability :
202 cases in stock

**** Non-US users please [sign in](#) or [get a quote](#) to view the proper price for your country. ****

Features:

- Suitable for long term tissue culture
- Manufactured in a class 100,000 clean room
- Frame made from virgin polystyrene.
- German high quality cover glass of superior optical quality, glass thickness is 0.170±0.005mm
- A USP class VI adhesive is used to assemble the cover glass and the plate.
- Sterilized by Gamma radiation.

Suitable for:

- Differential Interference Contrast (DIC)
- Widefield Fluorescence
- Confocal Microscopy
- Two-Photon and Multiphoton Microscopy
- Fluorescence Recovery After Photobleaching (FRAP)
- Förster Resonance Energy Transfer (FRET)
- Fluorescence Lifetime Imaging Microscopy (FLIM)
- Total Internal Reflection Fluorescence (TIRF)
- Super-Resolution Microscopy

Recommended for:

- Confocal Microscopy
- Super-Resolution Microscopy


Technical specifications

» [View technical specification of different coverslips.](#)

Coverslip	#1.5 high performance cover glass (0.170±0.005mm)
Length	57.00 mm

Width	25.00 mm
Height	10.80 mm
Temperature Range	-20°C to 50°C

Dimension diagram (units in mm)

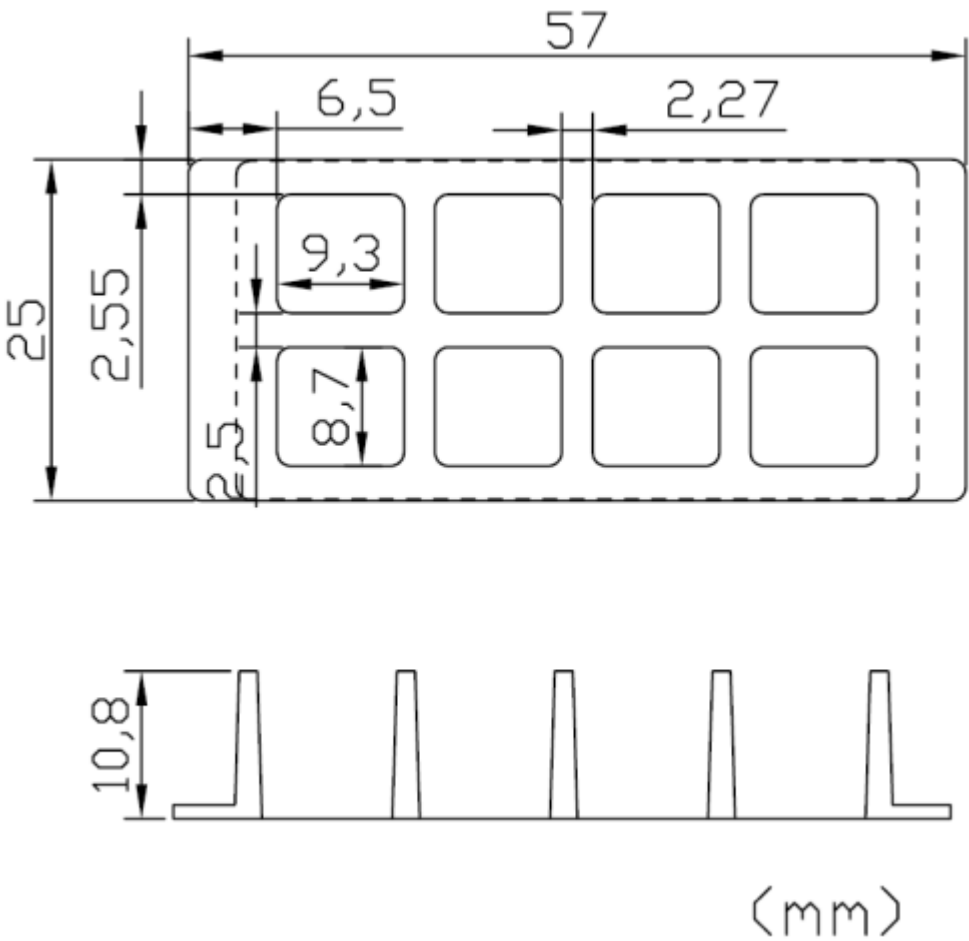

Latest cited publications on bioRxiv

- [TDP-43 pathology links innate and adaptive immunity in amyotrophic lateral sclerosis](#)
 Baggio A. Evangelista, et al., *bioRxiv - Cell Biology* 2024
 Quote: ... glass bottom slides (Cellvis, **C8-1.5H-N**), or glass-bottom culture dishes (Nunc ...
- [Synergistic block of SARS-CoV-2 infection by combined drug inhibition of the host entry factors PIKfyve kinase and TMPRSS2 protease](#)
 Alex J.B. Kreutzberger, et al., *bioRxiv - Microbiology* 2021
 Quote: ... 25 mm (Cellvis cat. **C8-1.5H-N**), Polydimethylsiloxane (PDMS ...
- [Mammary Epithelial Migration is EMT-Independent](#)
 Jing Chen, et al., *bioRxiv - Cell Biology* 2024
 Quote: ... or 8-well chamber (Cellvis, **C8-1.5H-N**), and beads and organoids were placed approximately 100 μm apart and covered with basic medium (DMEM/F12 or RPMI 1640 ...
- [The Arp2/3 complex is required for in situ haptotactic response of microglia to iC3b](#)
 Summer G. Paulson, et al., *bioRxiv - Neuroscience* 2025
 Quote: 8 well chamber dishes (Cellvis, **C8-1.5H-N**) were prepared with a coating of 10μg/mL fibronectin (Gibco ...
- [YTHDF1 mediates translational control by m6A mRNA methylation in adaptation to environmental challenges](#)
 Zhuoyue Shi, et al., *bioRxiv - Neuroscience* 2024
 Quote: 8 chambered cover glass systems (Cellvis **C8-1.5H-N**) were first prepared by coating them with 0.1 mg/mL poly-D-lysine (Sigma-Aldrich ...
- [Distinct transcriptomic profile of satellite cells contributes to preservation of neuromuscular junctions in extraocular muscles of ALS mice](#)
 Ang Li, et al., *bioRxiv - Cell Biology* 2024
 Quote: ... 8-well chambered cover glasses (Cellvis **C8-1.5H-N**) were coated with poly-L-ornithine and laminin as described above ...
- [Cathepsin B causes trogocytosis-mediated CAR T cell dysfunction](#)
 Kenneth A. Dietze, et al., *bioRxiv - Bioengineering* 2024
 Quote: Eight-well chambers (Cellvis, catalog no. **C8-1.5H-N**) were used for all experiments ...
- [RNA binding proteins and glycoRNAs form domains on the cell surface for cell penetrating peptide entry](#)
 Jonathan Perr, et al., *bioRxiv - Cell Biology* 2023
 Quote: ... in a 8-chambered coverglass system (Cellvis, **C8-1.5H-N**), treated with 500 μL fresh media ...
- [Transcription regulates bleb formation and stability independent of nuclear rigidity](#)
 Isabel K. Berg, et al., *bioRxiv - Cell Biology* 2022

Quote: ... Cells were plated in 8 well coverslips (Cellvis, **C8-1.5H-N**) grown and left untreated or treated with VPA and/or alpha amanitin ...

- [Temporal Dynamics and Stoichiometry in Notch Signaling - from Notch Synaptic Complex Formation to NICD Nuclear Entry](#)

Lena Tveriakhina, et al., *bioRxiv - Cell Biology* 2023

Quote: ... and seeded onto 8-well cover slips (Cellvis, **C8-1.5H-N**) in imaging media at 37°C in presence of 5% CO₂ at densities chosen to reach 30-50% confluency at the time of imaging the following day ...

[View all 63 references on labshake.com](#)

Cited Publications before 2019 (10)

- [The effects of IKK-beta inhibition on early NF-kappa-B activation and transcription of downstream genes](#)

MJ Bloom, et al., *Cellular Signalling* Volume 55, March 2019, Pages 17-25

Quote: "One day before the stimulation experiments, cells were seeded in 8-well chambered cover glass (Cellvis, Mountain View, CA) in 200 µl of culture medium at an estimated 35,000 cells/well."

- [The AP-1 complex regulates AXL expression and determines sensitivity to PI3Ka inhibition in esophagus and head and neck squamous cell carcinoma](#)

Mai Bdarny, et al., *bioRxiv*, October 08, 2018

Quote: "For IF, cells were seeded on 8-well glass slides (Cellvis, Cat no. C8-1.5H-N) for 48 hours"

- [Visualization of Calcium Ion Loss from Rotavirus during Cell Entry](#)

Eric N. Salgado, et al., *Journal of Virology*, 2018

Quote: "For live-cell imaging, BSC-1 cells were grown on 8-chambered number 1.5 cover glass slides (Cellvis) in a volume of 400 µl supplemented DMEM and allowed to grow overnight to 50% confluence."

- [EWI-2 Inhibits Cell-Cell Fusion at the Virological Synapse](#)

EE Whitaker, et al., *BioRxiv*, April 06, 2018

Quote: "Three days post infection, 3×10^5 infected cells were plated onto each well of 8-well glass-bottom plates (CellVis, Mountain View, CA, Cat. #C8-1.5HN) coated with 1:10 poly-L-Lysine in ddH₂O"

- [Quantitative Super-Resolution Imaging of Small RNAs in Bacterial Cells](#)

S Park, et al., *Bacterial Regulatory RNA* pp 199-212

Quote: "The imaging buffer: 10 mM NaCl, 50 mM Tris and 10% glucose in 1× PBS (for immunostained samples) or 4× SSC (for FISH samples) , with pH adjusted to 8. 5. 8-Well chambered cover glass (Nunc Lab-Tek 155409 or Cellvis C8-1.5HN)"

- [Analysis of microtubule dynamics heterogeneity in cell culture](#)

Anara Serikbaeva, et al., *Methods in Molecular Biology, Cellular Heterogeneity* pp 181-204

Quote: "For life imaging of transfected cells, we used eight-chambered cover glass system #1.5 high-performance cover glass (0.17 ± 0.005 mm) (In Vitro Scientific—www.invitrosoci.com)"

- [Size-dependent segregation controls macrophage phagocytosis of antibody-opsonized targets](#)

Matthew H Bakalar, et al., *BioRxiv*, January 19, 2018

Quote: "For imaging interfaces between cells and supported lipid bilayer coated beads, cells were seeded into 8-well imaging chambers with a cover-glass bottom (Cellvis) and beads were added to the wells once the cells had fully adhered to the cover-glass"

- [A Förster Resonance Energy Transfer Ratiometric Probe Based on Quantum Dot-Cresyl Violet for Imaging Hydrogen Sulfide in Living Cells](#)

BAI Min, et al., *Chinese Journal of Analytical Chemistry* Volume 46, Issue 1, January 2018, Pages 39-47

Quote: "Cell imaging analysis of H₂S in living cells using QDS-N3 ratiometric probe HeLa and MCF-7 cells were seeded in a 8-well Chambered Cover Glass (Cellvis, USA) at 40000 cells per well at 37 °C for overnight."

- [IKKa regulates human keratinocyte migration through surveillance of the redox environment](#)

Thomas S. Lissa, et al., *Journal of Cell Science*

Quote: "For RNA and scratch analyses, cells were cultured in 12-well plates (Corning, 3513) and 8-chamber glass-bottom dishes (In Vitro Scientific, C8-1.5H-N), respectively."

- [Effect of labeling density and time post labeling on quality of antibody-based super resolution microscopy images](#)

Amy M. Bittel, et al., *Proc. SPIE* 9331

Quote: "Immunofluorescence Labeling of Microtubules SKBR3 cells were grown in 8-well chambered coverglass slides (In Vitro Scientific) to ~75% confluency"

Cellvis (formerly In Vitro Scientific), P.O.Box 390959, Mountain View, CA 94039

Email: contact@cellvis.com, Phone(toll free): 1-866-203-7860

© 2007 - 2025 Cellvis (formerly In Vitro Scientific), all rights reserved

[Privacy](#) | [Contact Us](#)